



## Fabius® & More

Dräger offers a wide range of software and hardware options for:

- A Fabius® *plus* 8606800
- B Fabius® *plus XL* 8608555
- C Fabius® *Tiro* 8606000
- D Fabius® *GS premium* 8607000
- E Fabius® *MRI* 8607300

Not all products and accessories are available worldwide. Please contact your local Dräger sales representative for further information.



## CONTENT

<b>ACCESSORIES FOR FABIUS® BASIC UNIT</b>	<b>04</b>	<b>VENTILATION ACCESSORIES</b>	<b>23</b>
1. Conversion kits wall, ceiling and trolley variants	04	1. Disposable breathing circuits	23
2. Software upgrades	04	2. Reusable breathing circuits	30
3. Auxiliary fresh gas outlet	04	3. Breathing system filters and HMEs	34
4. Oxygen therapy	05	4. Catheter mounts and airway connectors	38
5. Storage and writing tray	06	5. Breathing masks	41
6. Additional workplace illumination	06	6. Set2go	47
7. Additional power outlets	07	7. Resuscitators	48
8. Additional mounting rails	07		
9. Miscellaneous accessories	08	<b>SODA LIME</b>	<b>49</b>
10. Accessories for IV	09	1. Disposable absorber	49
		2. Soda lime	49
<b>ACCESSORIES FOR BREATHING SYSTEM AND VENTILATOR</b>	<b>10</b>	<b>MONITORING COMBINATIONS</b>	<b>50</b>
1. Compact breathing system COSY	10		
2. Accessories for Compact breathing system cosy	10	<b>ENDOTRACHEAL ASPIRATION</b>	<b>52</b>
3. Upgrade Kit "Advanced low flow kit" for COSY	11	1. Ejector suction	52
4. Accessories for ventilator	12	2. Vacuum suction	53
5. Ventilator hoses	13	3. Accessories for endotracheal aspiration	54
6. Hoses for breathing system pneumatics	13		
7. Absorber and accessories	14		
8. Non-rebreathing systems	14		
<b>GAS MEASUREMENT</b>	<b>15</b>		
1. O <sub>2</sub> Measurement	15		
2. Multigas Measurement	15		
<b>CONSUMABLES FOR GAS AND FLOW MEASUREMENT</b>	<b>17</b>		
1. Gas Measurement	17		
2. Flow Measurement	17		
<b>VAPORIZERS</b>	<b>18</b>		
<b>MEDIA SUPPLY</b>	<b>19</b>		
1. Central gas supply	19		
2. Anesthetic gas scavenging systems (AGS)	20		

**ACCESSORIES FOR FABIUS® BASIC UNIT**

Dräger offers a wide range of software and hardware options for Fabius® *plus*, *plus XL*, *Tiro*, *GS premium* and *MRI*.

**ACCESSORIES FOR FABIUS® BASIC UNIT**

**1. CONVERSION KITS WALL, CEILING AND TROLLEY VARIANTS**

**Conversion kits wall, ceiling and trolley variants**

- Conversion Fabius <i>plus</i> : wall/ceiling mount to trolley	8608567
- Conversion Fabius <i>plus</i> : trolley to wall/ceiling mount	8608568
- Conversion Fabius <i>plus XL</i> : wall/ceiling mount to high trolley	2607081
- Conversion Fabius <i>plus XL</i> : wall/ceiling mount to trolley	2607082
- Conversion Fabius <i>plus XL</i> : trolley to wall/ceiling mount	2607083
- Conversion Fabius <i>Tiro</i> : wall/ceiling mount to trolley	8608570
- Conversion Fabius <i>Tiro</i> : trolley to wall/ceiling mount	8608569
- Conversion Fabius <i>GS premium</i> : trolley to ceiling mount	G19072

**2. SOFTWARE UPGRADES**

<b>A Software upgrades</b>	<b>8606410</b>
The conversion kit offers various software options and packages:	
- Pressure Controlled Ventilation (PCV)	
- Pressure Support (PS)	
- Synchronized Volume Controlled Ventilation with Pressure Support (SIMV/PS)	

**3. AUXILIARY FRESH GAS OUTLET**

There are different variants of external fresh gas outlets (A-cone) available for connecting non-rebreathing systems to Fabius®. These variants include left and right hand side versions with or without switch function. The switch versions allow users to simply switch between the Fabius® compact breathing system and the external fresh gas outlet.

<b>Upgrade kit auxiliary fresh gas outlet (A-cone) left</b>	<b>8607985</b>
<b>B Upgrade kit auxiliary fresh gas outlet (A-cone) right</b>	<b>8607986</b>
<b>C Upgrade kit auxiliary fresh gas outlet switch (A-cone switch) left</b>	<b>8607987</b>
<b>Upgrade kit auxiliary fresh gas outlet switch (A-Cone switch) right</b>	<b>8607988</b>

FABIUS® PLUS  
 FABIUS® PLUS XL  
 FABIUS® TIRO  
 FABIUS® GS PREMIUM  
 FABIUS® MRI



D-917-2010



MT-1341-2008



MT-3879-2008

## ACCESSORIES FOR FABIUS® BASIC UNIT

### 4. OXYGEN THERAPY

There are two options for connection an external O<sub>2</sub> therapy:

<b>A</b>	<b>Upgrade kit external O<sub>2</sub> flowtube</b> Additional oxygen flowmeter.	<b>8607694</b>	■	■	■
<b>B</b>	<b>Electronic O<sub>2</sub> flowmeter DigiFlow®, 16 L (541.02 fl. oz.), rail, NIST, 9/16" for sterile water</b> 16 L/min (541.02 fl. oz./min) oxygen flow meter with digital display of set flow. For directly plugging into an oxygen supply coupler. Adjustment from 0.3 to 3 L/min (10.14 to 101.44 fl. oz./min) in 0.1 increments, from 3 L/min (101.44 fl. oz./min) in 0.5 increments. With 9/16" port for sterile water systems or nozzle. Gas inlet: rotatable NIST union.	<b>MP01228</b>	■	■	■
<b>Mounting components for DigiFlow:</b>					
<b>C</b>	<b>Rail for GCX profile, 10 cm (3.94 in)</b>	<b>4118555</b>	■	■	■
	<b>NIST T-piece O<sub>2</sub></b>	<b>M36250</b>	■	■	■
	<b>Connecting hose O<sub>2</sub> NIST ISO</b>	<b>M36053</b>	■	■	■
<b>D</b>	<b>Adapter, with 6 mm (0.24 in) nozzle, e.g., for disposable breathing masks, 9/16" connection</b>	<b>MP00416</b>	■	■	■

FABIUS® PLUS  
FABIUS® PLUS XL  
FABIUS® TIRO  
FABIUS® GS PREMIUM  
FABIUS® MRI



## ACCESSORIES FOR FABIUS® BASIC UNIT

### 5. STORAGE AND WRITING TRAY



<b>A</b>	<b>Pivotable side tray</b> For mounting at the side rail. Load capacity max. 14 kg (30.86 lb).	<b>4117510</b>
<b>B</b>	<b>Pull-out writing tray</b> Additional pull-out writing tray. Load capacity max. 20 kg (44.09 lb).	<b>8608056</b>
<b>C</b>	<b>Flip-up tray</b> For mounting at the side rail. Load capacity max. 20 kg (44.09 lb).	<b>4118306</b>
	<b>Flip up tray for Fabius plus XL</b> For mounting at the side rail. Load capacity max. 20 kg (44.09 lb).	<b>2606709</b>
<b>D</b>	<b>Hook</b> Hook for accessories, e.g. breathing bag.	<b>M26349</b>
<b>E</b>	<b>Basket</b> For storage of accessories. Dimensions W × D × H: 30.5 × 20.3 × 15.2 cm (12 × 8 × 6 in).	<b>8606438</b>

### 6. ADDITIONAL WORKPLACE ILLUMINATION

With goseneck for flexible positioning.

<b>F</b>	<b>Halogen lamp, Euro</b>	<b>8604119</b>
<b>F</b>	<b>Halogen lamp, British Standard</b>	<b>8604120</b>
<b>F</b>	<b>Halogen lamp, US Standard</b>	<b>8604122</b>
	<b>Holder for halogen lamp</b> For mounting the halogen lamp to the side rail of the basic unit.	<b>8607019</b>

FABIUS® PLUS  
FABIUS® PLUS XL  
FABIUS® TIRO  
FABIUS® GS PREMIUM  
FABIUS® MRI

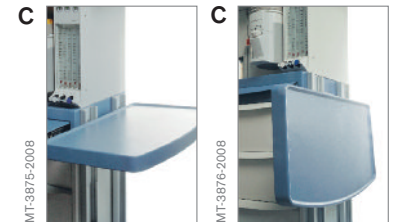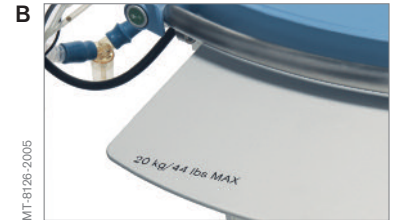

## ACCESSORIES FOR FABIUS® BASIC UNIT

### 7. ADDITIONAL POWER OUTLETS

<b>A</b>	<b>External power outlet strip</b>	<b>8607840</b>
	Four additional power sockets mounted on the mounting rail of Fabius.	



### 8. ADDITIONAL MOUNTING RAILS

Accessories can be fixed on the additional standard rails at the side or rear panel of Fabius® workstations.

	<b>Upgrade kit mounting rail, length 618 mm</b>	<b>8606259</b>
	For mounting at the rear panel of Fabius® GS <i>premium</i> and at the side panel of Fabius® Tiro, Fabius® GS <i>premium</i> .	



	<b>Upgrade kit mounting rail, length 567 mm</b>	<b>2605346</b>
	For mounting at the side panel of Fabius® <i>plus</i> . Alternatively the rail can be mounted at the rear panel with the Fabius® <i>plus</i> Pin Index version or in combination with a single cylinder holder.	



<b>B</b>	<b>Upgrade kit mounting rails, 2 pcs.</b>	<b>2606910</b>
	Rear panel, for mounting of gas cylinders or other accessories. Length: 340 mm (13.39 in).	

<b>C</b>	<b>Upgrade kit mounting rails – for control box, side front, 2 pcs.</b>	<b>2606918</b>
	Length: 266 mm (10.47 in).	



<b>D</b>	<b>Upgrade kit mounting rails – right side, low trolley, long/short, 2 pcs.</b>	<b>2606920</b>
	1 × long: length 1,186 mm (46.69 in) 1 × short: length 567 mm (22.32 in).	

	<b>Upgrade kit mounting rails – right side, high trolley, long/short, 2 pcs.</b>	<b>2606927</b>
	1 × long: length 1,246 mm (49.06 in) 1 × short: length 627 mm (24.69 in).	

## ACCESSORIES FOR FABIUS® BASIC UNIT

### 9. MISCELLANEOUS ACCESSORIES

<b>A</b>	<b>Earthing cable, 0.8 m (2.6 ft)</b>	<b>8301348</b>
	<b>Earthing cable, 3.2 m (10.5 ft)</b>	<b>8301349</b>
	<b>Earthing cable, 6 m (19.68 ft)</b>	<b>M27769</b>
	<b>Grips, 2 pcs.</b> Eases maneuverability of all Fabius® trolley variants. Two grips for mounting at the side rail.	<b>4117509</b>
	<b>Grips for Fabius® MRI, 2 pcs.</b>	<b>8607596</b>
<b>B</b>	<b>Hinged arm for hose and cable management</b> Helps to avoid cable clutter at the anesthesia workplace.	<b>M34080</b>
<b>C</b>	<b>Hinged arm holder</b> For mounting the hinged arm as well as the universal monitor arm.	<b>4117642</b>
	<b>CASTrGARD wheel guard device (large), set of 4</b>	<b>M36050</b>
	<b>CASTrGARD wheel guard device (small), set of 4</b> The CASTrGARDs are designed to prevent cables being overrun by the wheels of the device.	<b>M36049</b>
<b>D</b>	<b>Kit Cover mounting rails</b>	<b>4117441</b>
	<b>Adhesive Tape 40 mT/400 Gauss</b> To mark the floor of the MRI suite.	<b>8607593</b>

FABIUS® PLUS  
FABIUS® PLUS XL  
FABIUS® TIRO  
FABIUS® GS PREMIUM  
FABIUS® MRI





## ACCESSORIES FOR FABIUS® BASIC UNIT

### 10. ACCESSORIES FOR IV

<b>Kit pump and IV holder</b>	<b>8601863</b>
Two poles for mounting infusion bags and infusion pumps on the side rail.	
<b>A Pole for mounting IV pumps and hooks for infusion bags, height-adjustable</b>	<b>2606912</b>
<b>A Mounting arm for 2606912</b> Length: 300 mm (11.81 in).	<b>G25855</b>
<b>Down-pole for IV pumps</b>	<b>2606711</b>
<b>B Pole with hooks for infusion bags, non height adjustable</b>	<b>2606710</b>
<b>Mounting arm for 2606711 and 2606710</b> Length: 203 mm (7.99 in).	<b>2607097</b>

FABIUS® PLUS  
FABIUS® PLUS XL  
FABIUS® TIRO  
FABIUS® GS PREMIUM  
FABIUS® MRI



**ACCESSORIES FOR BREATHING SYSTEM AND VENTILATOR**

**ACCESSORIES FOR BREATHING SYSTEM AND VENTILATOR**

**1. COMPACT BREATHING SYSTEM COSY**

The Fabius® compact breathing system (COSY) is available with and without integrated warmer.

<b>A</b>	<b>COSY 2.6</b>	<b>8605797</b>	■ ■ ■ ■ ■
	<b>COSY 2.6 with warmer</b>	<b>8605796</b>	■ ■ ■ ■ ■
	Breathing system with integrated warmer (without power supply).		
	<b>Sightglass</b>	<b>M09230</b>	■ ■ ■ ■ ■
	For inspiration/expiration valve in COSY compact breathing system.		
	<b>Valve plate</b>	<b>M23225</b>	■ ■ ■ ■ ■
	Ceramic, valve disc.		

**2. ACCESSORIES FOR COMPACT BREATHING SYSTEM COSY**

<b>B</b>	<b>Hook for gas sample line</b>	<b>2606460</b>	■ ■ ■ ■ ■
	Placing the gas sample line in the hook for a tidy workplace (can not be combined with the Flow Sensor Guard and Cover Hood Cosy).		
<b>C</b>	<b>Airway pressure gauge</b>	<b>8607383</b>	■ ■ ■ ■ ■
<b>D</b>	<b>Cover Hood COSY, complete (can not be combined with the Flow Sensor Guard and hook for sample line)</b>	<b>8607111</b>	■ ■ ■ ■ ■
<b>E</b>	<b>Flow Sensor Guard (can not be combined with the Cover Hood COSY and hook for sample line)</b>	<b>8607218</b>	■ ■ ■ ■ ■
	<b>Pressure Port Guard</b>	<b>8608016</b>	■ ■ ■ ■ ■
<b>F</b>	<b>Flexible bag arm</b>	<b>8605088</b>	■ ■ ■ ■ ■
<b>G</b>	<b>Adjustable bag arm</b>	<b>4118380</b>	■ ■ ■ ■ ■

FABIUS® PLUS  
 FABIUS® PLUS XL  
 FABIUS® TIRO  
 FABIUS® GS PREMIUM  
 FABIUS® MRI



Not all articles are available worldwide.

**ACCESSORIES FOR BREATHING SYSTEM AND VENTILATOR**


**3. UPGRADE KIT "ADVANCED LOW FLOW KIT" FOR COSY**

The "Advanced low flow kit" for COSY was designed to prevent the condensation of water in the breathing circuit which typically occurs during long low flow cases. Therefore the kit includes:

- breathing system warmer
- power supply
- water trap for the ventilator hose

<b>A</b>	<b>Upgrade kit "Advanced low flow kit" for COSY 2.6</b>	<b>8605898</b>
----------	---	----------------

FABIUS® PLUS  
 FABIUS® PLUS XL  
 FABIUS® TIRO  
 FABIUS® GS PREMIUM  
 FABIUS® MRI

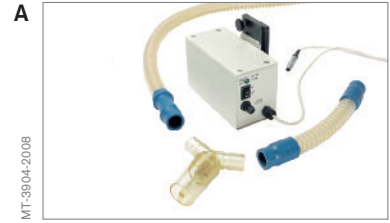


## ACCESSORIES FOR BREATHING SYSTEM AND VENTILATOR

### 4. ACCESSORIES FOR VENTILATOR

<b>A Ventilator diaphragm</b>	<b>2600650</b>
Ventilator diaphragm for Fabius® piston ventilator.	
<b>B Ventilator lid</b>	<b>8605828</b>

FABIUS® PLUS  
 FABIUS® PLUS XL  
 FABIUS® TIRO  
 FABIUS® GS PREMIUM  
 FABIUS® MRI

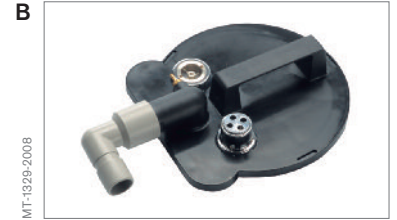

**ACCESSORIES FOR BREATHING SYSTEM AND VENTILATOR**

**5. VENTILATOR HOSES**

There are different lengths available for devices with left or right hand side COSY and for short or long COSY mounting arm.

<b>A</b>	<b>Ventilator hose, silicone, reusable, 1.1 m (3.6 ft), connector 22 mm/22 mm (0.87/0.87 in)</b>	<b>2166038</b>
<b>B</b>	<b>Ventilator hose, silicone, reusable, 1.5 m (4.9 ft), connector 22 mm/22 mm (0.87/0.87 in)</b>	<b>2166046</b>



**VENTILATOR HOSES FOR USE WITH WATERTRAP/LOW FLOW KIT**

<b>C</b>	<b>Ventilator hose, silicone, reusable, 0.6 m (23.6 in), connector 22 mm/22 mm (0.87/0.87 in)</b>	<b>2165627</b>
<b>D</b>	<b>Ventilator hose, silicone, reusable, 0.35 m (13.7 in), connector 22 mm/22 mm (0.87/0.87 in)</b>	<b>2165619</b>
	<b>Ventilator hose, silicone, reusable, 0.25 m (9.8 in), connector 22 mm/22 mm (0.87/0.87 in)</b>	<b>8604916</b>
<b>E</b>	<b>Watertrap</b> For use between the ventilator hoses.	<b>8404985</b>

**6. HOSES FOR BREATHING SYSTEM PNEUMATICS**

<b>F</b>	<b>Combination measurement hose PEEP-APL COSY left</b>	<b>8607055</b>
	<b>Combination measurement hose PEEP-APL COSY right</b>	<b>8607054</b>
	<b>Pressure measurement hose</b> Required length: 1.27 m (50 in) (COSY left) or 1.8 m (5.91 ft) (COSY right)	<b>1190520</b>
<b>G</b>	<b>Bacterial filter for pressure measurement hose, disposable</b> Bacterial filter to protect pressure sensor.	<b>8402868</b>

FABIUS® PLUS  
FABIUS® PLUS XL  
FABIUS® TIRO  
FABIUS® GS PREMIUM  
FABIUS® MRI

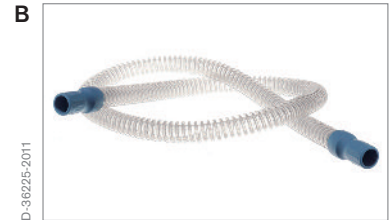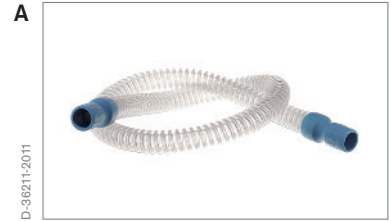

Not all articles are available worldwide.

## ACCESSORIES FOR BREATHING SYSTEM AND VENTILATOR

### 7. ABSORBER AND ACCESSORIES

<b>A</b>	<b>Reusable absorber container</b>	<b>M33719</b>
	<b>Absorber insert for reusable absorber</b>	<b>M33720</b>
<b>B</b>	<b>Absorber adapter CLIC, for Dräger anesthesia workstations</b>	<b>MX50090</b>
	To use the disposable absorber at the anesthesia workstations Fabius® Family.	

### 8. NON-REBREATHING SYSTEMS

<b>C</b>	<b>Medec Waters Set</b>	<b>M36037</b>
	Semi-open breathing system with integrated adjustable pressure valve and AGS connection.	
	<b>Adapter Waters Set AGS</b>	<b>M36086</b>
	<b>Catheter mount, 15 cm (5.91 in) (Waters Set)</b>	<b>MX21710</b>

FABIUS® PLUS  
 FABIUS® PLUS XL  
 FABIUS® TIRO  
 FABIUS® GS PREMIUM  
 FABIUS® MRI



## GAS MEASUREMENT

### GAS MEASUREMENT

#### 1. O<sub>2</sub> MEASUREMENT

<b>A</b>	<b>O<sub>2</sub> sensor housing</b>	<b>8606055</b>	
	For the integrated O <sub>2</sub> measurement of Fabius®. For O <sub>2</sub> cell see page 17.		

#### 2. MULTIGAS MEASUREMENT

	<b>Vamos</b>	<b>6871990</b>	
	Standalone gas monitor for the measurement of volatile agents, N <sub>2</sub> O and CO <sub>2</sub> . Variants: – CO <sub>2</sub> and SpO <sub>2</sub> – Anesthetic agents, CO <sub>2</sub> and N <sub>2</sub> O – Anesthetic agents, CO <sub>2</sub> , N <sub>2</sub> O and SpO <sub>2</sub> . Option: Vamos can be configured with rechargeable batteries.		

<b>B</b>	<b>Vamos Plus with Auto ID</b>	<b>6871990</b>	
	Vamos with automatic identification of the anesthetic agent and dual agent display when two agents are used. Variants: – Anesthetic agents, CO <sub>2</sub> and N <sub>2</sub> O – Anesthetic agents, CO <sub>2</sub> , N <sub>2</sub> O and SpO <sub>2</sub> . Option: Vamos can be configured with rechargeable batteries.		

#### Mounting options for Vamos

	There are different mounting options available in the Vamos configurator:		
	– Top of the anesthesia device		
	– Side arm		
	– Wall mounted rail system		

FABIUS® PLUS  
 FABIUS® PLUS XL  
 FABIUS® TIRO  
 FABIUS® GS PREMIUM  
 FABIUS® MRI



**GAS AND FLOW MEASUREMENT**

FABIUS® PLUS  
 FABIUS® PLUS XL  
 FABIUS® TIRO  
 FABIUS® GS PREMIUM  
 FABIUS® MRI

<p><b>A SCIO Four</b> <span style="float: right;"><b>6871810</b></span></p> <p>The Patient Gas Module SCIO Four is the anesthetic gas measurement module for Dräger patient monitors and available in different variants:</p> <ul style="list-style-type: none"> <li>- SCIO Four: Anesthetic agents, CO<sub>2</sub> and N<sub>2</sub>O</li> <li>- SCIO Four plus: Anesthetic agents, CO<sub>2</sub>, N<sub>2</sub>O, automatic identification of the anesthetic agent and mixture detection</li> <li>- SCIO Four Oxi: Anesthetic agents, CO<sub>2</sub>, N<sub>2</sub>O and O<sub>2</sub> (free of consumption)</li> <li>- SCIO Four Oxi plus: Anesthetic agents, CO<sub>2</sub>, N<sub>2</sub>O, automatic identification of the anesthetic agent, mixture detection and O<sub>2</sub> (free of consumption).</li> </ul>	
<p><b>Data Line Fabius®</b> <span style="float: right;"><b>8601473</b></span></p> <p>Connection cable for automatic Desflurane compensation of the flow measurement in Fabius®.</p>	
<p><b>Integration options for SCIO</b></p> <p>The Upgrade Kits allow the integration of the SCIO.</p>	
<p><b>B Fabius GS premium SCIO integration</b> <span style="float: right;"><b>8606137</b></span></p> <p>SCIO and power supplies are integrated into the left side of the trolley.</p>	
<p><b>C Fabius Tiro SCIO integration</b> <span style="float: right;"><b>2603900</b></span></p> <p>SCIO and power supplies are integrated instead of the upper drawer of the Fabius® Tiro trolley.</p>	
<p><b>D Fabius plus XL SCIO integration</b> <span style="float: right;"><b>2606923</b></span></p> <p>SCIO and power supplies are integrated instead of the upper drawer of the Fabius® plus XL trolley.</p>	
<p><b>Additional RS232 port</b> <span style="float: right;"><b>2606944</b></span></p> <p>A second RS232 port for simultaneous data export and automatic Desflurane compensation.</p>	





**GAS MEASUREMENT / FLOW SENSORS**

**CONSUMABLES FOR GAS AND FLOW MEASUREMENT**

**1. GAS MEASUREMENT**

<b>A</b>	<b>Watertrap WaterLock® 2, 12 pcs.</b> Replacement for watertrap 6870567.	<b>6872130</b>	■	■	■	■
<b>B</b>	<b>Adapter, set of 6 pcs.</b> To connect the new WaterLock® 2 for all Dräger devices with the holder for the old watertrap 6870267.	<b>6872128</b>	■	■	■	■
<b>C</b>	<b>Sampleline set, 2.5 m (8.2 ft), disposable, 10 pcs.</b>	<b>8290286</b>	■	■	■	■
	<b>Sampleline, 5 m (16.4 ft), 10 pcs.</b>	<b>8290602</b>	■	■	■	■
<b>D</b>	<b>Sample gas return kit for Fabius®</b> For recirculation of gas into the breathing system.	<b>M32692</b>	■	■	■	■

<b>E</b>	<b>Oxygen sensor, O<sub>2</sub> capsule</b> Mainstream oxygen sensor. Sensors are electrochemical systems intended for the measurement of O <sub>2</sub> partial pressure. The sensors work according to the galvanic cell principle. The gas mixture to be monitored diffuses through a plastic membrane into the sensor's liquid electrolyte. The electrolyte contains a working electrode (gold electrode) and a counter-electrode (lead electrode). The oxygen to be monitored is electrochemically reduced at the working electrode.	<b>6850645</b>	■	■	■	■
----------	--	----------------	---	---	---	---

**2. FLOW MEASUREMENT**

<b>F</b>	<b>Flow sensor Spirolog®, polycarbonate, disinfectable, 5 pcs.</b> Spirolog® Flowsensor is a hot-wire sensor for measuring volumetric gas flow (flow) delivered by Dräger anesthesia machines and ventilators. Spirolog® is not suitable for sterilization, but disinfection. Material: Acrylonitrile-Butadiene-Styrene (ABS).	<b>8403735</b>	■	■	■	■
----------	---	----------------	---	---	---	---

<b>G</b>	<b>Flow sensor SpiroLife®, PPSU, disinfectable and sterilizable</b> The reusable sensor of Dräger Medical. SpiroLife® is the only device which combines the advantages of the hot-wire anemometer technology with the possibility of sterilization. Material: Polycarbonate (PC).	<b>MK01900</b>	■	■	■	■
----------	--	----------------	---	---	---	---

FABIUS® PLUS  
FABIUS® PLUS XL  
FABIUS® TIRO  
FABIUS® GS PREMIUM  
FABIUS® MRI



Not all articles are available worldwide.

## VAPORIZERS

### VAPORIZERS

#### VAPORE 2000

<b>A</b>	<b>D-Vapor</b>	<b>M35500</b>	■	■	■	■
<b>B</b>	<b>Vapor 2000</b>	<b>M35054</b>	■	■	■	■
	Isoflurane					
	Sevoflurane					
	Halothane					

#### DRÄGER FILL™ FILLING ADAPTERS

	<b>Filling adapter, (i) isoflurane (violet)</b>	<b>M36110</b>	■	■	■	■
<b>C</b>	<b>Filling adapter, (s) sevoflurane (yellow)</b>	<b>M36120</b>	■	■	■	■
	<b>Filling adapter, (h) halothane (red)</b>	<b>M36090</b>	■	■	■	■

#### VAPORIZER STANDBY HOLDER

For mounting two vaporizers on the side rail.

<b>D</b>	<b>Vaporizer standby holder Dräger</b>	<b>4115837</b>	■	■	■	■
	<b>Vaporizer standby holder Selectatec</b>	<b>8607178</b>	■	■	■	■

FABIUS® PLUS  
 FABIUS® PLUS XL  
 FABIUS® TIRO  
 FABIUS® GS PREMIUM  
 FABIUS® MRI



**MEDIA SUPPLY**

**MEDIA SUPPLY**

**1. CENTRAL GAS SUPPLY**

**A Dräger supply hoses, selectable by the adjacent configuration number 8601697**

Design, color coding and connection thread of all Dräger supply hoses comply with the respective legal regulations in the indicated countries.


**ADAPTER HOSES**

**O<sub>2</sub> adapter hose for Australia M34446**

**N<sub>2</sub>O adapter hose for Australia M34447**

**AIR adapter hose for Australia M34448**

FABIUS® PLUS  
FABIUS® PLUS XL  
FABIUS® TIRO  
FABIUS® GS PREMIUM  
FABIUS® MRI



**A**

D-2286-2014



Not all articles are available worldwide.

**MEDIA SUPPLY**

**2. ANESTHETIC GAS RECEIVING SYSTEM (AGS)**

<b>A</b>	<b>Anesthetic Gas Scavenging System AGS</b>	<b>M33300</b>
	For scavenging excess gas from anesthetic workstations. For reducing the anaesthetic agent concentration in the vicinity of anaesthetic workstations. Operating range of suction: 30 to 50 L/min (7.93 to 13.21 Imp. gal).	
<b>B</b>	<b>NPT regulation valve</b>	<b>8605566</b>
	Only for USA/Canada/Costa Rica/Puerto Rico.	
	<b>Hose coupling</b>	<b>M33149</b>
	For USA/Canada.	
<b>C</b>	<b>Transfer hose adapter</b>	<b>M29430</b>
	Cone adapter 30 mm for mounting the transfer hose.	

FABIUS® PLUS  
 FABIUS® TIRO  
 FABIUS® GS  
 FABIUS® GS PREMIUM  
 FABIUS® MRI



**MEDIA SUPPLY**

**TRANSFER HOSE**

For connecting the Fabius® breathing system to the Anesthetic Gas Scavenging System (AGS).


<b>A</b>	<b>AGS transfer hose, 1 m (3.3 ft)</b>	<b>M33295</b>	■	■	■	■
<b>A</b>	<b>AGS transfer hose, 1.5 m (4.9 ft)</b>	<b>M33296</b>	■	■	■	■
	For Fabius® with breathing system on the right hand side and Fabius® MRI.					

**SUCTION HOSES AND CONNECTORS**

Between anesthetic gas scavenging system (AGS) and central gas supply.

<b>B</b>	<b>AGS suction hose, 1.5 m (4.9 ft)</b>	<b>M33297</b>	■	■	■	■
<b>B</b>	<b>AGS suction hose, 3 m (9.8 ft)</b>	<b>M33298</b>	■	■	■	■
<b>B</b>	<b>AGS suction hose, 5 m (16.4 ft)</b>	<b>M33299</b>	■	■	■	■
<b>C</b>	<b>AGS suction hose, British Standard, 3.5 m (11 ft)</b>	<b>M36051</b>	■	■	■	■
	Only for Great Britain/Ireland.					
<b>D</b>	<b>AGS scavenging hose connector Japan/Australia/New Zealand</b>	<b>M36048</b>	■	■	■	■

FABIUS® PLUS  
FABIUS® PLUS XL  
FABIUS® TIRO  
FABIUS® GS PREMIUM  
FABIUS® MRI

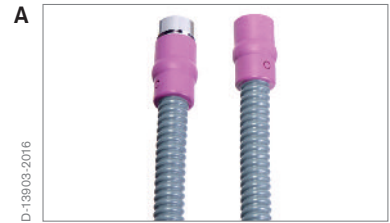

Not all articles are available worldwide.

**MEDIA SUPPLY**



**CONNECTOR FOR AGS**

<b>A</b>	<b>Anesthetic Gas Scavenging connector, 45° angled, DIN</b>	<b>G60440</b>
<b>B</b>	<b>Anesthetic Gas Scavenging connector, 45° angled, EN-ISO</b>	<b>G60495</b>
<b>C</b>	<b>Anesthetic Gas Scavenging connector, straight, EN-ISO</b>	<b>G60580</b>
<b>D</b>	<b>AGS Ejector</b> For use with AGS (M33300) if a central Anesthetic Gas Scavenging System is not available. The ejector generates a vacuum for an active scavenging of excess anesthetic and expiration gas from the AGS. A fresh gas hose for connecting a T-piece is included.	<b>M36175</b>
<b>E</b>	<b>Adapter 3-Way AIR NIST</b> 3-way splitter for handling the AGS ejector.	<b>M36056</b>

**MOUNTING OF AGS**

<b>F</b>	<b>Hook for AGS</b> For mounting the AGS at the standard side rails of the trolley.	<b>4117515</b>
	<b>AGS rail claw</b> For mounting the AGS at the wall rail.	<b>M32967</b>

FABIUS® PLUS  
FABIUS® PLUS XL  
FABIUS® TIRO  
FABIUS® GS PREMIUM  
FABIUS® MRI

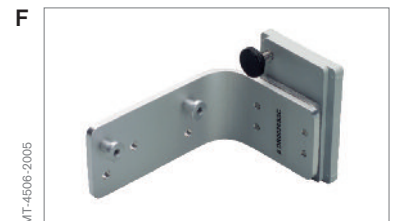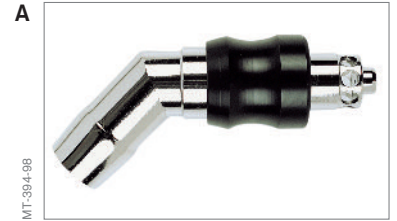

**VENTILATION ACCESSORIES**



**VENTILATION ACCESSORIES**

**1. DISPOSABLE BREATHING CIRCUITS**

**DISPOSABLE BREATHING CIRCUITS FOR ADULTS**

<b>A</b>	<b>Breathing circuit Vent Set, latex-free, disposable, basic, 1.5 m (4.92 ft), 25 pcs.</b>	<b>MP00300</b>
<b>B</b>	<b>Breathing circuit VentStar®, latex-free, disposable, basic, 1.8 m (5.91 ft), 25 pcs.</b>	<b>MP00349</b>
<b>C</b>	<b>Breathing circuit VentStar®, latex-free, disposable, basic, 2.5 m (8.2 ft), 10 pcs.</b>	<b>MP00350</b>
<b>D</b>	<b>Breathing circuit Vent Set, latex-free, disposable, extendable, 1.8 m (5.91 ft), 25 pcs.</b>	<b>MP00305</b>
<b>E</b>	<b>Breathing circuit VentStar®, latex-free, disposable, extendable, 2.2 m (7.22 ft), 25 pcs.</b>	<b>MP00355</b>

FABIUS® PLUS  
 FABIUS® PLUS XL  
 FABIUS® TIRO  
 FABIUS® GS PREMIUM  
 FABIUS® MRI

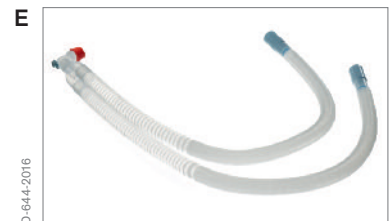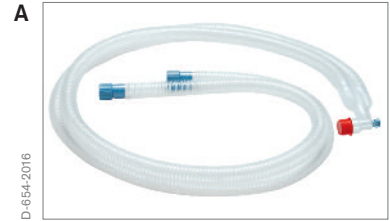


Not all articles are available worldwide.

**VENTILATION ACCESSORIES**

FABIUS® PLUS  
 FABIUS® PLUS XL  
 FABIUS® TIRO  
 FABIUS® GS PREMIUM  
 FABIUS® MRI

<b>A</b>	<b>Breathing circuit Vent Set, latex-free, disposable, coaxial, 1.5 m (4.92 ft), 20 pcs.</b>	<b>MP00315</b>
<b>B</b>	<b>Breathing circuit Vent Set, latex-free, disposable, coax, with holder, 1.5 m (4.92 ft) 20 pcs.</b>	<b>MP00316</b>
<b>C</b>	<b>Breathing circuit VentStar®, latex-free, disposable, water trap, 1.8 m (5.91 ft), 25 pcs.</b>	<b>MP00337</b>
<b>D</b>	<b>Breathing circuit VentStar®, latex-free, disposable, smoothbore, for MRI, 3 m (9.8 ft) , 10 pcs.</b>	<b>MP00327</b>

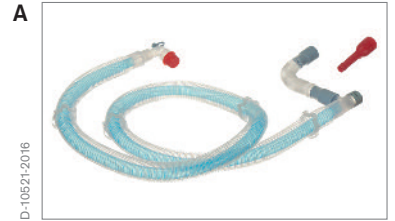




**VENTILATION ACCESSORIES**

FABIUS® PLUS  
FABIUS® PLUS XL  
FABIUS® TIRO  
FABIUS® GS PREMIUM  
FABIUS® MRI


<b>A</b>	<b>Breathing bag set VentStar®, latex-free, disposable, 2 L (67.63 fl. oz.), 0.8 m (2.6 ft), 25 pcs.</b>	<b>MP00222</b>
<b>B</b>	<b>Breathing bag set VentStar®, latex-free, disposable, 2 L (67.63 fl. oz.), 1.2 m (3.9 ft), 25 pcs.</b>	<b>MP00380</b>
<b>B</b>	<b>Breathing bag set VentStar®, latex-free, disposable, 2 L (67.63 fl. oz.), 1.5 m (4.92 ft), 25 pcs.</b>	<b>MP00381</b>
<b>C</b>	<b>Anesthesia set, latex-free, disposable, basic, 1.5 m/0.8 m (4.92/2.6 ft), 25 pcs.</b>	<b>MP00301</b>
<b>D</b>	<b>Anesthesia set, latex-free, disposable, extendable, 1.8 m/0.8 m (5.91/2.6 ft), 25 pcs.</b>	<b>MP00303</b>
<b>E</b>	<b>Anesthesia set VentStar®, latex-free, disposable, w/o Luer Lock, 1.8 m/0.8 m (5.91/2.6 ft), 25 pcs.</b>	<b>MP00330</b>

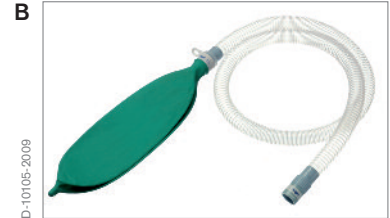

FABIUS® PLUS  
 FABIUS® PLUS XL  
 FABIUS® TIRO  
 FABIUS® GS PREMIUM  
 FABIUS® MRI

## VENTILATION ACCESSORIES

<b>A</b> Anesthesia set VentStar®, latex-free, disposable, basic, 2 L (67.63 fl. oz.), 1.8 m/1.5 m (5.91/4.92 ft), 25 pcs.	<b>MP00365</b>
<b>B</b> Anesthesia set VentStar®, latex-free, disposable, basic, 2 L (67.63 fl. oz.), 2.5 m/1.5 m (8.2/4.92 ft), 10 pcs.	<b>MP00366</b>
<b>C</b> Anesthesia set VentStar®, latex-free, disposable, WT, w/o LL, 2 L (67.63 fl. oz.), 1.8 m/1.5 m (5.91/4.92 ft), 10 pcs.	<b>MP00371</b>
<b>D</b> Anesthesia set VentStar®, latex-free, disposable, WT, 2 L (67.63 fl. oz.), 1.8 m/1.5 m (5.91/4.92 ft), 10 pcs.	<b>MP00372</b>


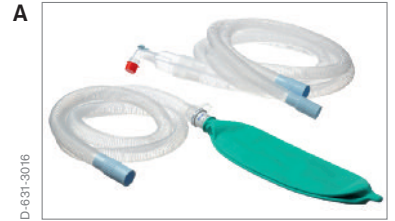


**VENTILATION ACCESSORIES**

**DISPOSABLE BREATHING CIRCUITS FOR PEDIATRICS**

<b>A</b>	<b>Breathing circuit VentStar®, latex-free, disposable, P, basic, w/o LL, 1.5 m (4.92 ft), 25 pcs.</b>	<b>MP00338</b>
<b>B</b>	<b>Breathing circuit VentStar®, latex-free, disposable, P, basic, 1.8 m (5.91 ft), 25 pcs.</b>	<b>MP00351</b>
<b>C</b>	<b>Breathing circuit VentStar®, latex-free, disposable, P, basic, 2.5 m (8.2 ft), 25 pcs.</b>	<b>MP00352</b>
<b>D</b>	<b>Breathing circuit VentStar®, latex-free, disposable, P, water trap, 1.8 m (5.91 ft), 25 pcs.</b>	<b>MP00361</b>
<b>E</b>	<b>Breathing circuit VentStar®, latex-free, disposable, P, coaxial, 1.5 m (4.92 ft), 20 pcs.</b>	<b>MP00379</b>

FABIUS® PLUS  
 FABIUS® PLUS XL  
 FABIUS® TIRO  
 FABIUS® GS PREMIUM  
 FABIUS® MRI

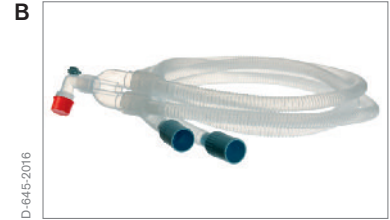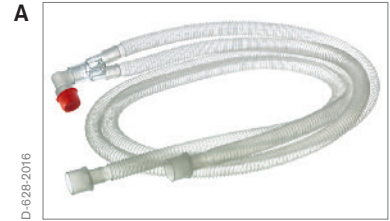

Not all articles are available worldwide.

FABIUS® PLUS  
 FABIUS® PLUS XL  
 FABIUS® TIRO  
 FABIUS® GS PREMIUM  
 FABIUS® MRI

## VENTILATION ACCESSORIES

<b>A Breathing bag set VentStar®, latex-free, disposable, P, 1 L (33.81 fl. oz.), 1.1 m (3.61 ft), 25 pcs.</b>	<b>MP00384</b>
<b>B Anesthesia set VentStar®, latex-free, disposable, P, w/o Luer Lock, 1.5 m/1.1 m (4.92/3.61 ft), 25 pcs.</b>	<b>MP00331</b>
<b>C Anesthesia set VentStar®, latex-free, disposable, P, WT, w/o LL, 1 L (33.81 fl. oz.), 1.8 m/1.1 m (5.91/3.61 ft), 10 pcs.</b>	<b>MP00375</b>
<b>D Anesthesia set VentStar®, latex-free, disposable, P, WT, 1 L (33.81 fl. oz.), 1.8 m/1.1 m (5.91/3.61 ft), 10 pcs.</b>	<b>MP00374</b>



**VENTILATION ACCESSORIES**

**DISPOSABLE BREATHING CIRCUITS FOR NEONATES**

<b>A</b>	<b>Breathing circuit VentStar®, latex-free, disposable, N, MRI, 3 m (9.8 ft) , 10 pcs.</b>	<b>MP00329</b>
<b>B</b>	<b>Breathing bag set VentStar®, latex-free, disposable, N, 0.5 L (16,91 fl. oz.), 1.1 m (3.61 ft), 25 pcs.</b>	<b>MP00383</b>
<b>C</b>	<b>Anesthesia set VentStar®, latex-free, disposable, N, 1.8 m/1.1 m (5.91/3.61 ft), 25 pcs.</b>	<b>MP00333</b>

FABIUS® PLUS  
 FABIUS® PLUS XL  
 FABIUS® TIRO  
 FABIUS® GS PREMIUM  
 FABIUS® MRI




Not all articles are available worldwide.

**VENTILATION ACCESSORIES**

**2. REUSABLE BREATHING CIRCUITS**

**REUSABLE BREATHING CIRCUITS FOR ADULTS AND PEDIATRICS**

<b>A</b>	<b>BlueSet Anesthesia (P/A)</b>	<b>M30146</b>
Consists of:		
①	Breathing hose, silicone, reusable, 1.5 m (4.92 ft), connector 22 mm/22 mm (0.87/0.87 in)	<b>2166046</b>
②	Breathing hose, silicone, reusable, 1.1 m (3.61 ft), connector 22 mm/22 mm (0.87/0.87 in)	<b>2166038</b>
③	Breathing bag, silicone, reusable, 2.3 L (77,77 fl. oz.)	<b>2166062</b>
④	Y-piece, polyamide, reusable, A, with angled Luer Lock	<b>M33278</b>
⑤	ISO double connector, reusable, 22 mm/22 mm (0.87/0.87 in)	<b>M25647</b>

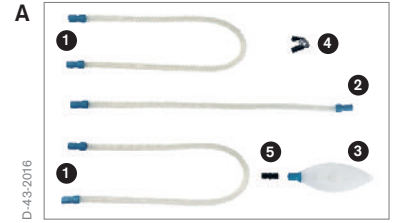


D-43 2016

**VENTILATION ACCESSORIES**

**REUSABLE BREATHING CIRCUITS FOR NEONATES**

<b>A BlueSet Anesthesia (N)</b>	<b>M27542</b>
Consists of:	
<b>1 Breathing hose, silicone, reusable, P, 1.5 m (4.92 ft), connector 10 mm/22 mm (0.39/0.87 in), Ø 10 mm (0.39 in)</b>	<b>2165961</b>
<b>2 Breathing hose, silicone, reusable, P, 1.1 m (3.61 ft), connector 22 mm/22 mm (0.87/0.87 in), Ø 10 mm (0.39 in)</b>	<b>2165678</b>
<b>3 Breathing bag, silicone, reusable, 0.5 L (16,91 fl. oz.)</b>	<b>2166054</b>
<b>4 ISO double connector, reusable, 22 mm/22 mm (0.87/0.87 in)</b>	<b>M25647</b>
<b>5 Y-piece, reusable, P, angled, with suction opening</b>	<b>8403075</b>
<b>6 Luer Lock elbow, N</b>	<b>M34568</b>
<b>6 Cap Luer Lock</b>	<b>M34461</b>
<b>7 Cap set, reusable, 7 mm (0,28 in), 5 pcs.</b>	<b>8402953</b>

FABIUS® PLUS  
 FABIUS® PLUS XL  
 FABIUS® TIRO  
 FABIUS® GS PREMIUM  
 FABIUS® MRI

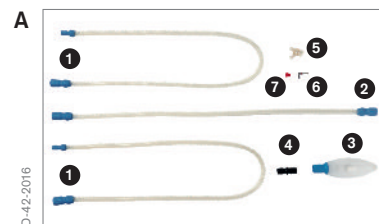

Not all articles are available worldwide.

## VENTILATION ACCESSORIES

### Y-PIECES AND CONNECTORS

<b>A</b>	<b>Y-piece, polyamide, reusable, A, straight</b>	<b>M25650</b>
<b>B</b>	<b>Y-piece, reusable, P, straight</b>	<b>M27077</b>
<b>C</b>	<b>Y-piece, polyamide, reusable, A, angled, without Luer Lock</b>	<b>M25682</b>
<b>D</b>	<b>Y-piece, polyamide, reusable, A, with angled Luer Lock</b>	<b>M33278</b>
<b>E</b>	<b>Y-piece, reusable, P, angled, with suction opening</b>	<b>8403075</b>
<b>F</b>	<b>ISO double connector, reusable, 22 mm/22 mm (0.87/0.87 in)</b>	<b>M25647</b>

FABIUS® PLUS  
 FABIUS® PLUS XL  
 FABIUS® TIRO  
 FABIUS® GS PREMIUM  
 FABIUS® MRI






**VENTILATION ACCESSORIES**

**BREATHING BAGS, REUSABLE**

<b>A</b>	Transparent breathing bags of silicone elastomer, latex free				
<b>1</b>	<b>Breathing bag, silicone, reusable, 0.5 L (17.5 fl. oz.)</b>	<b>2166054</b>	■	■	■
<b>2</b>	<b>Breathing bag, silicone, reusable, 1.5 L (52.8 fl. oz.)</b>	<b>2166127</b>	■	■	■
<b>3</b>	<b>Breathing bag, silicone, reusable, 2.3 L (77,77 fl. oz.)</b>	<b>2166062</b>	■	■	■
<b>4</b>	<b>Breathing bag, silicone, reusable, 3 L (101.44 fl. oz.)</b>	<b>2166070</b>	■	■	■

FABIUS® PLUS  
 FABIUS® PLUS XL  
 FABIUS® TIRO  
 FABIUS® GS PREMIUM  
 FABIUS® MRI

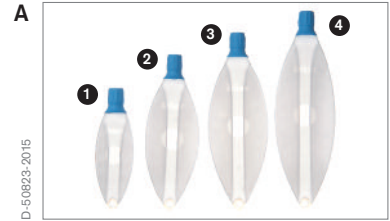

Not all articles are available worldwide.

**VENTILATION ACCESSORIES**

**3. BREATHING SYSTEM FILTERS AND HMEs**

**ELECTROSTATIC FILTERS AND HMEs**

The TwinStar® filter/HME from Dräger is a combination of an efficient HME and a high-performance breathing system filter, thereby significantly contributing to infection prophylaxis in ventilation treatment.

- Combination of efficient HME and high-performance breathing system filter
- High retention rates for bacteria and viruses
- High return of heat and moisture
- Reduced condensation in the breathing circuit
- Minimal workload and low cost
- Sampling connector with tethered cap
- Sampling port in convenient 45° angle
- Transparent housing for visual control
- Standardized connectors for safe connection to other components
- Clear labeling and blue color coding for quick identification

<b>A</b>	<b>Electrostatic filter and HME TwinStar® 90, disposable, 90 mL (3.04 fl. oz.) deadspace, 50 pcs.</b>	<b>MP01800</b>	■	■	■	■
<b>B</b>	<b>Electrostatic filter and HME TwinStar® 65A, disposable, 65 mL (2.2 fl. oz.) deadspace, angled, 50 pcs.</b>	<b>MP01810</b>	■	■	■	■
<b>C</b>	<b>Electrostatic filter and HME TwinStar® 55, disposable, 55 mL (1.86 fl. oz.) deadspace, 50 pcs.</b>	<b>MP01805</b>	■	■	■	■
<b>D</b>	<b>Electrostatic filter and HME TwinStar® 25, disposable, 25 mL (0.85 fl. oz.) deadspace, 50 pcs.</b>	<b>MP01815</b>	■	■	■	■
<b>E</b>	<b>Electrostatic filter and HME TwinStar® 10A, disposable, 10 mL (0.34 fl. oz.) deadspace, angled, 50 pcs.</b>	<b>MP01825</b>	■	■	■	■
<b>F</b>	<b>Electrostatic filter and HME TwinStar® 8, disposable, 8 mL (0.27 fl. oz.) deadspace, 50 pcs.</b>	<b>MP01820</b>	■	■	■	■
<b>G</b>	<b>Mechanical filter and HME TwinStar® HEPA, disposable, 55 mL (1.86 fl. oz.) deadspace, 50 pcs.</b>	<b>MP01801</b>	■	■	■	■

FABIUS® PLUS  
 FABIUS® PLUS XL  
 FABIUS® TIRO  
 FABIUS® GS PREMIUM  
 FABIUS® MRI



Not all articles are available worldwide.

**VENTILATION ACCESSORIES**

**ELECTROSTATIC FILTERS**

The CareStar® breathing system filter significantly contributes to infection prophylaxis in anesthesia and ventilation treatment through high retention rates.

- High-performance electrostatic filter
- Patient-side or device-side use
- Minimal workload and low cost
- Sampling connector with tethered cap
- Sampling port in convenient 45° angle
- Transparent housing for visual control
- Standardized connectors for safe connection to other components
- Clear labeling and red color coding for quick identification

<b>A</b>	<b>Electrostatic filter CareStar® 45, disposable, 45 mL (1.52 fl. oz.) deadspace, 50 pcs.</b>	<b>MP01755</b>				
<b>B</b>	<b>Electrostatic filter CareStar® 40A, disposable, 40 mL (1.35 fl. oz.) deadspace, angled, 50 pcs.</b>	<b>MP01765</b>				
<b>C</b>	<b>Electrostatic filter CareStar® 30, disposable, 30 mL (1.01 fl. oz.) deadspace, 50 pcs.</b>	<b>MP01770</b>				

FABIUS® PLUS  
 FABIUS® PLUS XL  
 FABIUS® TIRO  
 FABIUS® GS PREMIUM  
 FABIUS® MRI



Not all articles are available worldwide.

## VENTILATION ACCESSORIES

### MECHANICAL FILTERS

The SafeStar® mechanical high performance breathing system filter shows excellent retention rates, thereby significantly contributing to infection prophylaxis in anesthesia and ventilation treatment.

- High-performance mechanical HEPA filter
- High hydrophobicity
- Patient-side or device-side use
- Sampling connector with tethered cap
- Sampling port in convenient 45° angle
- Transparent housing for visual control
- Standardized connectors for safe connection to other components
- Clear labeling and red color coding for quick identification

<b>A</b>	<b>Mechanical filter SafeStar® 80, disposable, 80 mL (2.71 fl. oz.) deadspace, 50 pcs.</b>	<b>MP01785</b>
<b>B</b>	<b>Mechanical filter SafeStar® 60A, disposable, 60 mL (2.03 fl. oz.) deadspace, angled, 50 pcs.</b>	<b>MP01795</b>
<b>C</b>	<b>Mechanical filter SafeStar® 55, disposable, 55 mL (1.86 fl. oz.) deadspace, 50 pcs.</b>	<b>MP01790</b>

FABIUS® PLUS  
 FABIUS® PLUS XL  
 FABIUS® TIRO  
 FABIUS® GS PREMIUM  
 FABIUS® MRI

**A**

MT-4386-2007



**B**

MT-4746-2007



**C**

MT-1165-2006



**VENTILATION ACCESSORIES**

**HMEs**

HumidStar® heat and moisture exchangers from Dräger partially take over the function of the upper airways and efficiently heat and humidify the inspired air.

- High return of heat and moisture
- Reduced condensation in the breathing circuit
- Minimal workload and low cost
- Sampling connector with tethered cap
- Sampling port in convenient 45° angle
- Transparent housing for visual control
- Standardized connectors for safe connection to other components
- Clear labeling and green color coding for quick identification

<b>A</b>	<b>HME HumidStar® 55, disposable, 55 mL (1.86 fl. oz.) deadspace, 50 pcs.</b>	<b>MP01730</b>	■	■	■
<b>B</b>	<b>HME HumidStar® 25, disposable, 25 mL (0.85 fl. oz.) deadspace, 50 pcs.</b>	<b>MP01735</b>	■	■	■
<b>C</b>	<b>HME HumidStar® 10A, disposable, 10 mL (0.34 fl. oz.) deadspace, angled, 50 pcs.</b>	<b>MP01740</b>	■	■	■
<b>D</b>	<b>HME HumidStar® 2, disposable, 2 mL (0,07 fl. oz.) deadspace, 50 pcs.</b>	<b>MP01745</b>	■	■	■

FABIUS® PLUS  
FABIUS® PLUS XL  
FABIUS® TIRO  
FABIUS® GS PREMIUM  
FABIUS® MRI



Not all articles are available worldwide.

## VENTILATION ACCESSORIES

### 4. CATHETER MOUNTS AND AIRWAY CONNECTORS

- Smoothbore products provide low resistance to flow while external coil reduces kinks and blockages
- Extendable product provides a lightweight and flexible solution
- Double swivel connector increases ease of use
- Double cap to maintain PEEP during suctioning or bronchoscopy
- Each catheter mount and airway connector is equipped with a red safety cap so that the system is completely closed and protected until directly prior to use on the patient
- 22F/15M connector on device side as a two-in-one solution
- Standardized connectors for safe connection to other components
- Luer Lock port for gas sampling
- Disposable convenience enhances workflow
- PVC-free design of gas leading components eliminates risk of softeners such as DEHP

<b>A</b>	<b>Catheter mount ErgoStar® CM 40, disposable, 10 cm (3.94 in), double swivel elbow, 50 pcs.</b>	<b>MP01840</b>
<b>B</b>	<b>Catheter mount ErgoStar® CM 45, disposable, 15 cm (5.91 in), double swivel elbow, 50 pcs.</b>	<b>MP01845</b>
<b>C</b>	<b>Catheter mount ErgoStar® CM 50, disposable, 10 – 16 cm (3.94 – 6.3 in), extendable, 50 pcs.</b>	<b>MP01850</b>
<b>D</b>	<b>Catheter mount ErgoStar® CM 55, disposable, 15 cm (5.91 in), straight, 50 pcs.</b>	<b>MP01855</b>
<b>E</b>	<b>Catheter mount ErgoStar® CM 60, disposable, 10 cm (3.94 in), small, double swivel elbow, 50 pcs.</b>	<b>MP01860</b>
<b>F</b>	<b>Airway connector ErgoStar® AC 90, disposable, double-swivel elbow w/ double cap, 100 pcs.</b>	<b>MP01890</b>



FABIUS® PLUS  
FABIUS® PLUS XL  
FABIUS® TIRO  
FABIUS® GS PREMIUM  
FABIUS® MRI

## VENTILATION ACCESSORIES

### COMBINATIONS OF FILTER AND TUBE EXTENTION

<b>A</b>	<b>CombiStar Filter HME flex</b>	<b>MP04240</b>
	Combination of flexible tube extension MP01850 and filter TwinStar® 55 (MP01805).	
	<b>CombiStar F-HME HEPA flex</b>	<b>MP04242</b>
	Combination of flexible tube extension MP01850 and filter TwinStar® HEPA (MP01801).	
<b>B</b>	<b>CombiStar Filter flex</b>	<b>MP04244</b>
	Combination of flexible tube extension MP01850 and filter CareStar® 30 (MP01770).	
<b>C</b>	<b>CombiStar mech. Filter flex</b>	<b>MP04246</b>
	Combination of flexible tube extension MP01850 and filter SafeStar® 55 (MP01790).	



## VENTILATION ACCESSORIES

<b>A</b>	<b>CombiStar Filter HME straight</b> Combination of straight tube extention MP01855 and filter TwinStar® 55 (MP01805).	<b>MP04230</b>
	<b>CombiStar F-HME HEPA straight</b> Combination of straight tube extention MP01855 and filter TwinStar® HEPA (MP01801).	<b>MP04232</b>
<b>B</b>	<b>CombiStar Filter straight</b> Combination of straight tube extention MP01855 and filter CareStar® 30 (MP01770).	<b>MP04234</b>
<b>C</b>	<b>CombiStar mech. Filter straight</b> Combination of straight tube extention MP01855 and filter SafeStar® 55 (MP01790).	<b>MP04236</b>
<b>D</b>	<b>CombiStar A Filter HME straight</b> Combination of airway connector MP01845 und Filter TwinStar® 55 (MP01805).	<b>MP01671</b>







**VENTILATION ACCESSORIES**

**5. BREATHING MASKS**

**BREATHING MASKS, DISPOSABLE**

A	① Anesthesia mask LiteStar®, disposable, w/o hook ring, size 1, premat., 30 pcs.	MP01501	■	■	■
	② Anesthesia mask LiteStar®, disposable, w/o hook ring, size 2, baby, 30 pcs.	MP01502	■	■	■
	③ Anesthesia mask LiteStar®, disposable, w/o hook ring, size 3, infant, 30 pcs.	MP01503	■	■	■
	④ Anesthesia mask LiteStar®, disposable, w/o hook ring, size 4, adult S, 30 pcs.	MP01504	■	■	■
	⑤ Anesthesia mask LiteStar®, disposable, w/o hook ring, size 5, adult M, 30 pcs.	MP01505	■	■	■
	⑥ Anesthesia mask LiteStar®, disposable, w/o hook ring, size 6, adult L, 30 pcs.	MP01506	■	■	■
B	① Hook ring for LiteStar® anesthesia masks, disposable, 15 mm (0.59 in) OD, for sizes 1 – 2, 30 pcs.	MP01551	■	■	■
	② Hook ring for LiteStar® anesthesia masks, disposable, 22 mm (0.87 in) ID, for sizes 3 – 6, 30 pcs.	MP01552	■	■	■


FABIUS® PLUS  
 FABIUS® PLUS XL  
 FABIUS® TIRO  
 FABIUS® GS PREMIUM  
 FABIUS® MRI

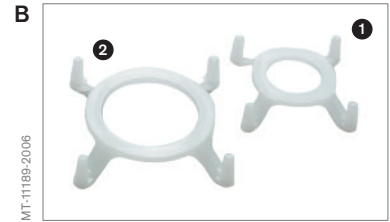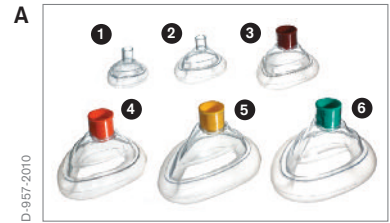

**VENTILATION ACCESSORIES**

A	<b>1</b> Anesthesia mask ComfortStar®, disposable, hook rg, std., size 0, premat., 20 pcs.	MP01510
	<b>2</b> Anesthesia mask ComfortStar®, disposable, hook rg, std., size 1, neo., 20 pcs.	MP01511
	<b>3</b> Anesthesia mask ComfortStar®, disposable, hook rg, std., size 2, baby, 20 pcs.	MP01512
	<b>4</b> Anesthesia mask ComfortStar®, disposable, hook rg, std., size 3, inf., 20 pcs.	MP01513
	<b>5</b> Anesthesia mask ComfortStar®, disposable, hook rg, std., size 4, adt. S, 20 pcs.	MP01514
	<b>6</b> Anesthesia mask ComfortStar®, disposable, hook rg, std., size 5, adt. M, 20 pcs.	MP01515
	<b>7</b> Anesthesia mask ComfortStar®, disposable, hook rg, std., size 6, adt. L, 20 pcs.	MP01516



FABIUS® PLUS  
 FABIUS® PLUS XL  
 FABIUS® TIRO  
 FABIUS® GS PREMIUM  
 FABIUS® MRI

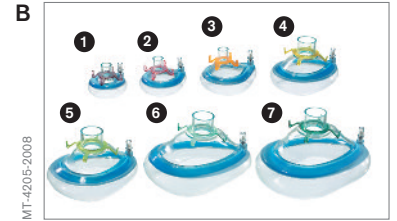

FABIUS® PLUS  
FABIUS® PLUS XL  
FABIUS® TIRO  
FABIUS® GS PREMIUM  
FABIUS® MRI

**VENTILATION ACCESSORIES**




A	① Anesthesia mask ComfortStar®, disposable, hook rg, strawb., size 1, neo., 20 pcs.	MP01521	■	■	■
	② Anesthesia mask ComfortStar®, disposable, hook rg, strawb., size 2, baby, 20 pcs.	MP01522	■	■	■
	③ Anesthesia mask ComfortStar®, disposable, hook rg, strawb., size 3, inf., 20 pcs.	MP01523	■	■	■
	④ Anesthesia mask ComfortStar®, disposable, hook rg, strawb., size 4, adt. S, 20 pcs.	MP01524	■	■	■
B	① Anesthesia mask ComfortStar®, disposable, hook rg, bubblg., size 1, neo., 20 pcs.	MP01531	■	■	■
	② Anesthesia mask ComfortStar®, disposable, hook rg, bubblg., size 2, baby, 20 pcs.	MP01532	■	■	■
	③ Anesthesia mask ComfortStar®, disposable, hook rg, bubblg., size 3, inf., 20 pcs.	MP01533	■	■	■
	④ Anesthesia mask ComfortStar®, disposable, hook rg, bubblg., size 4, adt. S, 20 pcs.	MP01534	■	■	■
C	① Anesthesia mask ComfortStar®, disposable, hook rg, mint, size 3, inf., 20 pcs.	MP01543	■	■	■
	② Anesthesia mask ComfortStar®, disposable, hook rg, mint, size 4, adult S, 20 pcs.	MP01544	■	■	■
	③ Anesthesia mask ComfortStar®, disposable, hook rg, mint, size 5, adult M, 20 pcs.	MP01545	■	■	■
	④ Anesthesia mask ComfortStar®, disposable, hook rg, mint, size 6, adult L, 20 pcs.	MP01546	■	■	■

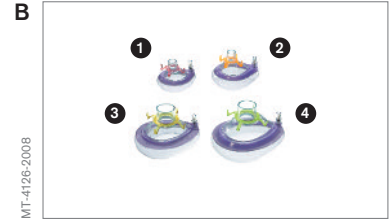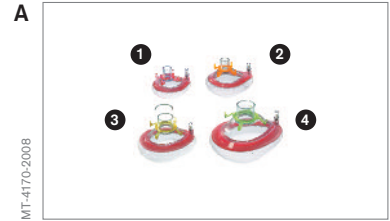

## VENTILATION ACCESSORIES

A	① Headgear, reusable, silicone, small	MP01557
	② Headgear, reusable, silicone, large	MP01558
B	① Headgear, disposable, soft fabric, small	MP01559
	② Headgear, disposable, soft fabric, large	MP01560
C	① Headgear, disposable, soft rubber, small	MP01561
	② Headgear, disposable, soft rubber, large	MP01562
D	Two-sided pump ball for air management	MP01590



FABIUS® PLUS  
 FABIUS® PLUS XL  
 FABIUS® TIRO  
 FABIUS® GS PREMIUM  
 FABIUS® MRI

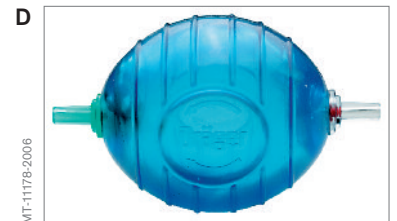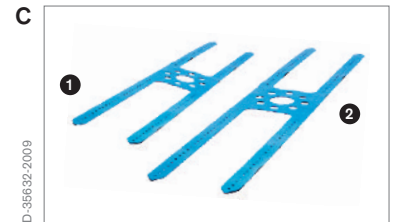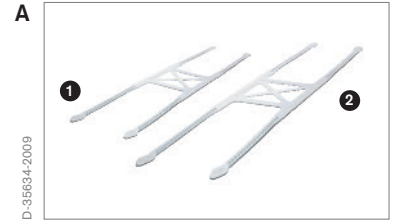

**VENTILATION ACCESSORIES**

**BREATHING MASKS, REUSABLE**

A	① Silicone mask, reusable, size 00, premature baby	MP02900
	② Silicone mask, reusable, size 0, small neonate	MP01910
	③ Silicone mask, reusable, size 1, large neonate	MP01911
	④ Silicone mask, reusable, size 2, small child	MP01912
	⑤ Silicone mask, reusable, size 3, large child	MP01913
	⑥ Silicone mask, reusable, size 4, small adult	MP01914
	⑦ Silicone mask, reusable, size 5, large adult	MP01915
B	① Hook ring for silicone masks, reusable, 15 mm (0.59 in) OD, for sizes 0 – 1, 10 pcs.	MP01916
	② Hook ring for silicone masks, reusable, 22 mm (0.87 in) ID, for sizes 2 – 5, 10 pcs.	MP01917

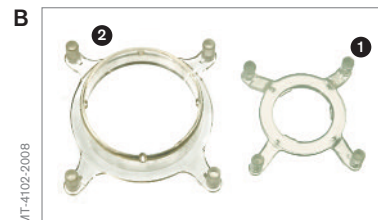
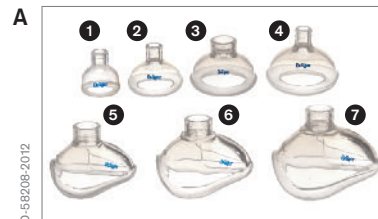


**VENTILATION ACCESSORIES**

**OXYGEN AND AEROSOL THERAPY MASKS**

<p><b>A</b> O2Star® Non-rebreather oxygen mask with connection tube 2.1 m (6.89 ft), disposable, S, 50 pcs.</p>	<p><b>MP01930</b></p>
<p><b>B</b> O2Star® Non-rebreather oxygen mask with connection tube 2.1 m (6.89 ft), disposable, M/L, 50 pcs.</p>	<p><b>MP01920</b></p>
<p><b>C</b> O2Star® Medium concentration oxygen mask with connection tube 2.1 m (6.89 ft), disposable, S, 50 pcs.</p>	<p><b>MP01931</b></p>
<p><b>D</b> O2Star® Medium concentration oxygen mask with connection tube 2.1 m (6.89 ft), disposable, M/L, 50 pcs.</p>	<p><b>MP01921</b></p>

FABIUS® PLUS  
 FABIUS® PLUS XL  
 FABIUS® TIRO  
 FABIUS® GS PREMIUM  
 FABIUS® MRI

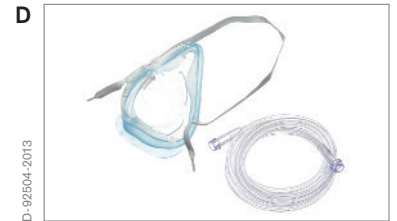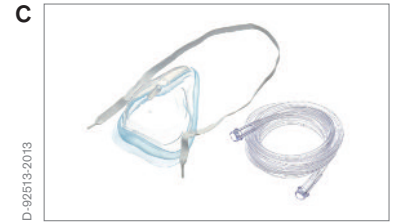



## VENTILATION ACCESSORIES

### 7. RESUSCITATORS

<b>A</b>	<b>Man. breathing bag MR-100, reusable, latex-free, A, w/ pressure valve, 60±10 mbar</b>	<b>5702881</b>	■	■	■
<b>B</b>	<b>Man. breathing bag MR-100, reusable, latex-free, P, w/ pressure valve, relief &gt;40±5 mbar</b>	<b>5702321</b>	■	■	■
<b>C</b>	<b>Man. breathing bag MR-100, reusable, latex-free, N, w/ pressure valve, relief &gt;40±5 mbar</b>	<b>5702322</b>	■	■	■
<b>D</b>	<b>Hook</b> Hook for accessories, e.g. breathing bag.	<b>M26349</b>	■	■	■
<b>E</b>	<b>O<sub>2</sub> Oxidem 3000 demand valve</b> For administration of O <sub>2</sub> during ventilation with a manual ventilation bag.	<b>2M86296</b>	■	■	■

FABIUS® PLUS  
FABIUS® PLUS XL  
FABIUS® TIRO  
FABIUS® GS PREMIUM  
FABIUS® MRI





**SODA LIME**

**SODA LIME**

**1. DISPOSABLE ABSORBER**

<b>A</b>	<b>CO<sub>2</sub> absorber CLIC 800+, disposable absorber, 1.2 L (40.58 fl. oz.), 6 pcs.</b>	<b>MX00004</b>
<b>B</b>	<b>CO<sub>2</sub> absorber CLIC Free, disposable absorber, 1.2 L (40.58 fl. oz.), 6 pcs.</b>	<b>MX50100</b>
<b>C</b>	<b>Absorber adapter CLIC, for Dräger anesthesia workstations</b>	<b>MX50090</b>

**2. SODA LIME**

Soda lime in pellet form with color indicator from white to violet. For the absorption of carbon dioxide (CO<sub>2</sub>) from breathing gas in inhalation anesthesia machines with rebreathing system.

Note: Fabius® devices are only authorised for use with Dräger CO<sub>2</sub> absorbers.

<b>D</b>	<b>Soda lime Drägersorb® 800+, for CO<sub>2</sub> absorption, canister, 5 L (169.07 fl. oz.), 2 pcs.</b>	<b>MX00001</b>
<b>E</b>	<b>Soda lime Drägersorb® Free, for CO<sub>2</sub> absorption, canister, 5 L (169.07 fl. oz.), 2 pcs.</b>	<b>MX50050</b>
<b>F</b>	<b>Canister opener, for canisters MX00001 and MX50050, 5 pcs.</b>	<b>MX10209</b>

FABIUS® PLUS  
FABIUS® PLUS XL  
FABIUS® TIRO  
FABIUS® GS PREMIUM  
FABIUS® MRI



## MONITORING COMBINATIONS

### MONITORING COMBINATIONS


All Fabius® devices have open interface architecture and can be combined with hemodynamic monitoring. In combination with Dräger Infinity® patient monitors and the SCIO patient gas module the Fabius® Workstation provides integrated gas monitoring as well as enhanced ventilation monitoring including spirometry. Dräger offers CE-certified combinations for many hemodynamic monitors.

Infinity® monitors and SCIO can not be used with the Fabius® MRI. Cooperations with third party monitoring suppliers exist for the Fabius MRI (see next page).


Ask your sales representative for an updated list of approved mounting and monitoring combinations. He may use the Acombi tool for visualisation.



D-11244-2016

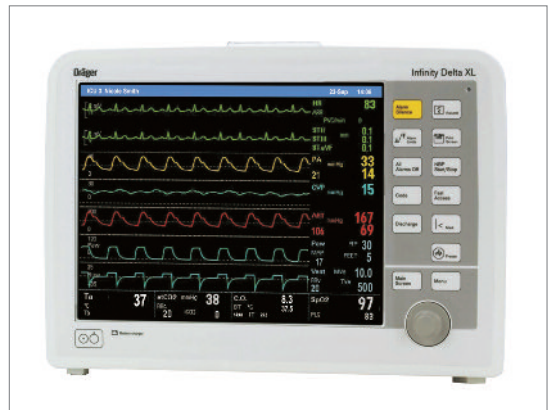

D-36662-2009

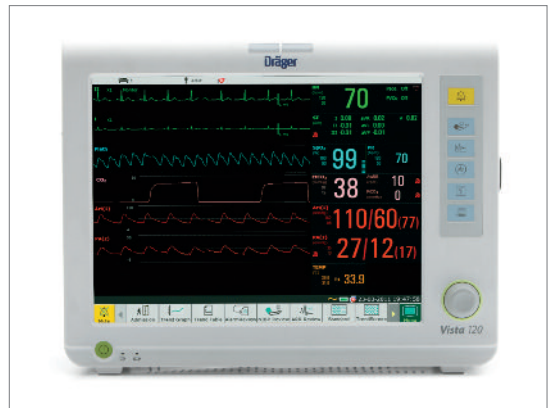

D-6829-2014



D-12436-2016

## MONITORING COMBINATIONS

### FABIUS MRI WITH TESLA M3

Fabius MRI can be combined with MIPM Tesla M3 patient monitor for hemodynamic monitoring in the MRI suite. A wireless remote monitor displays ventilation parameter and alarms from Fabius MRI in the control room.

Ask your sales representative for details.



D-48891-2015



D-48896-2015

**ENDOTRACHEAL ASPIRATION**

**ENDOTRACHEAL ASPIRATION**

The endotracheal suction is used to extract secretion from the airway, e.g. for ventilated patients.

**1. EJECTOR SUCTION**

For connection to central supply system or to 3-way AIR/O<sub>2</sub> adaptor on the gasinlet of Fabius®. Not for US.

**A Ejector suction NIST for bronchial suction, gas connection AIR** **MK01420**

**A Ejector suction NIST AIR O<sub>2</sub> for bronchial suction, gas connection AIR/O<sub>2</sub>** **MK01418**

**SUCTION ACCESSORIES**

**B Adapter 3-Way AIR NIST** **M36056**  
For the gas supply of the suction ejector.

**Adapter 3-Way NIST AIR/O<sub>2</sub>** **4117391**

**Hose NIST AIR, 1 m (3.3 ft) (connection to device)** **M34948**

**C Hose NIST AIR/O<sub>2</sub>, 1.5 m (4.9 ft)** **M33145**

**D Ejector VarioAir® B compact, for bronchial suction, gas connection AIR** **MP00536**  
Including jar set.

**Holder for portable ejector** **8608029**

FABIUS® PLUS  
FABIUS® PLUS XL  
FABIUS® TIRO  
FABIUS® GS PREMIUM  
FABIUS® MRI



## ENDOTRACHEAL ASPIRATION

### 2. VACUUM SUCTION

**A Vacuum suction NIST for bronchial suction, gas connection VAC** **MK01422**

**A Vacuum suction DISS for bronchial suction, gas connection VAC** **MK03140**

**A Vacuum suction DISS Canada for bronchial suction, gas connection VAC** **MK03320**

### MOUNTING

**B Rail for GCX profile, 10 cm (3.94 in)** **4118555**  
 For mounting the ejector aspirator at the standard rail of Fabius®.

FABIUS® PLUS  
 FABIUS® PLUS XL  
 FABIUS® TIRO  
 FABIUS® GS PREMIUM  
 FABIUS® MRI



## ENDOTRACHEAL ASPIRATION

### 3. ACCESSORIES FOR ENDOTRACHEAL ASPIRATION

<b>A</b>	<b>Secretion jar set with clamp</b> Including suction accessories.	<b>MP00646</b>
<b>B</b>	<b>Swivel arm for secretion jar set</b> For mounting the secretion jar at the standard rail of Fabius®.	<b>4118453</b>
	<b>MRI suction Kit GCX</b> Holder for certain non Dräger secretion jars.	<b>8608100</b>


FABIUS® PLUS  
FABIUS® PLUS XL  
FABIUS® TIRO  
FABIUS® GS PREMIUM  
FABIUS® MRI

**A**



D-44720-2015

**B**

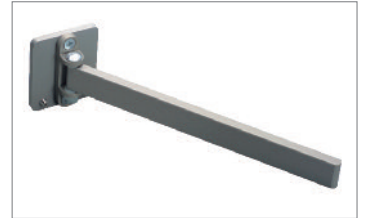

MT-4540-2005

**ENDOTRACHEAL ASPIRATION**

FABIUS® PLUS  
 FABIUS® PLUS XL  
 FABIUS® TIRO  
 FABIUS® GS PREMIUM  
 FABIUS® MRI

**A Bacterial filter, for endotr. suction for Primus®, Perseus®, Fabius®, Zeus® IE, disposable, 20 pcs. MK00514**

You can use the bacterial filter to reliably protect the suction devices at anesthesia workplaces against contamination.

Technical data:  
 Bacterial retention: 99.999%\*\*  
 Viral retention: 99.999%\*\*  
 Material: Acryl  
 Recommended duration of use: 7 days  
 Filtration method: Hydrophobic (HEPA)\*\*

**B Bacterial and viral filter system VarioSafe, for VarioAir®, VarioVac®, VarioO2, disposable, 10 pcs. MP00555**

The new VarioSafe® disposable filter system reliably protects the VarioAir®, VarioVac® and VarioO2 series and the patient environment against contamination.

Technical data:  
 Maximum capacity: 90 mL (3.04 fl. oz.)  
 Bacterial retention: 99.99998%\*  
 Viral retention: 99.99998%\*  
 Material: PS (polystyrol)  
 Recommended duration of use: 7 days  
 Filtration method: Mechanical and hydrophobic (HEPA)\*

\*Test has been performed by Nelson Lab, USA, before 2010.

\*\*According to EN 1822-1 Filter for suspended particles (HEPA and ULPA).




## ENDOTRACHEAL ASPIRATION

### DISPOSABLE SUCTION ACCESSORIES

<p><b>A Secretion jar VacuSmart®, PE (Polyethylene), disposable, 700 mL (23.67 fl. oz.), 25 pcs.</b></p>	<p><b>MX23023</b></p>
<p>The VacuSmart® is a disposable cartridge and has been specially designed and thoroughly tested to for compatibility with the Dräger endotracheal suction unit using 700 mL (23.67 fl. oz.) secretion jars, so there is no need to modify the suction unit. It is adaptable within few seconds. An overflow safeguard reliably protects the aspirator and the environment against contamination.</p>	
<p>Technical data:                  Maximum capacity: 700 mL (23.67 fl. oz.)                  Maximum capacity up to the overflow safeguard: 1,000 mL (33,81 fl. oz.)                  Overflow safeguard: Integrated, hydrophobic filter                  Recommended duration of use: 24 h</p>	
<p><b>B Secretion jar VacuSmart® Gel, PE (Polyethylene), disposable, 700 mL (23.67 fl. oz.), 25 pcs.</b></p>	<p><b>MX23025</b></p>
<p>The VacuSmart® Gel is a disposable cartridge with an integrated gelling agent for all Dräger endotracheal suction units using 700 mL (23.67 fl. oz.) secretion jars. The gelling agent coagulates the bronchial secretion and prevent thereby the leakage e.g. in case of disposal in a waste press. The integrated sodium polyacrylate is additionally packed in a water soluble bag in order to avoid skin contact.</p>	
<p>Technical data:                  Maximum capacity: 700 mL (23.67 fl. oz.)                  Maximum capacity up to the overflow safeguard: 1,000 mL (33,81 fl. oz.)                  Gelling agent: Sodium polyacrylate                  Overflow safeguard: Integrated, hydrophobic filter                  Recommended duration of use: 24 h</p>	
<p><b>C Suction tube VacuSmart® Tube, Polyethylene, disposable, 2 m, w/ finger tip, 25 pcs.</b></p>	<p><b>MX23026</b></p>
<p><b>D Suction set, disposable, VacuSmart®, 700 mL (23.67 fl. oz.) and VacuSmart® Tube, 25 pcs. each</b></p>	<p><b>MP00061</b></p>
<p><b>D Set VacuSmart® Gel and Tube, 25 pcs.</b></p>	<p><b>MP03211</b></p>

FABIUS® PLUS  
 FABIUS® PLUS XL  
 FABIUS® TIRO  
 FABIUS® GS PREMIUM  
 FABIUS® MRI

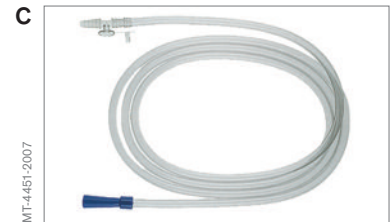






**ENDOTRACHEAL ASPIRATION**

**REUSABLE SUCTION ACCESSORIES**

FABIUS® PLUS  
 FABIUS® PLUS XL  
 FABIUS® TIRO  
 FABIUS® GS PREMIUM  
 FABIUS® MRI

<b>Reusable suction accessory set</b>	<b>8601179</b>	■	■	■	■
Consists of:					
<b>A Secretion jar, PSU (Polysulfone), reusable, 700 mL (23.67 fl. oz.), autoclavable up to 134°C</b>	<b>2M85594</b>	■	■	■	■
<b>B Secretion sight glass, 5 pcs.</b>	<b>M07582</b>	■	■	■	■
<b>C Cap, for 700 mL (23.67 fl. oz.) secretion jar, with overflow safeguard</b>	<b>2M85011</b>	■	■	■	■
<b>D Suction hose, silicone, reusable, 1.5 m (4.92 ft), for medical suction</b>	<b>M25780</b>	■	■	■	■
<b>E Connection hose, silicone, reusable, ID 5 mm (0.2 in), wall thickness 2 mm (0.08 in), by meter</b>	<b>1203606</b>	■	■	■	■
<b>F Elbow, angled hose connection for secretion jar caps</b>	<b>2M19063</b>	■	■	■	■

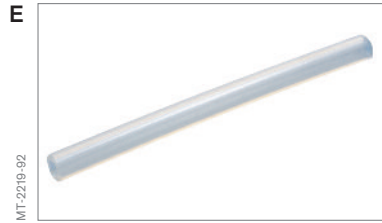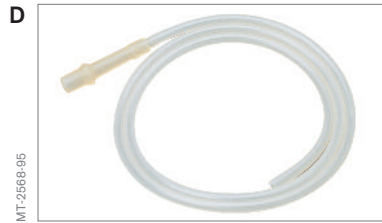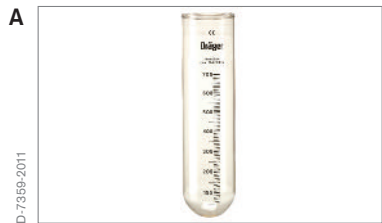

Not all articles are available worldwide.





Not all products, features, or services are for sale in all countries.  
Mentioned Trademarks are only registered in certain countries and not necessarily in the country in which this material is released. Go to [www.draeger.com/trademarks](http://www.draeger.com/trademarks) to find the current status.

#### CORPORATE HEADQUARTERS

Drägerwerk AG & Co. KGaA  
Moislinger Allee 53–55  
23558 Lübeck, Germany

[www.draeger.com](http://www.draeger.com)

#### Manufacturer:

Drägerwerk AG & Co. KGaA  
Moislinger Allee 53–55  
23558 Lübeck, Germany

Locate your Regional Sales  
Representative at:  
[www.draeger.com/contact](http://www.draeger.com/contact)



#### REGION DACH

Drägerwerk AG & Co. KGaA  
Moislinger Allee 53–55  
23558 Lübeck, Germany  
Tel +49 451 882 0  
Fax +49 451 882 2080  
[info@draeger.com](mailto:info@draeger.com)

#### REGION EUROPE

Drägerwerk AG & Co. KGaA  
Moislinger Allee 53–55  
23558 Lübeck, Germany  
Tel +49 451 882 0  
Fax +49 451 882 2080  
[info@draeger.com](mailto:info@draeger.com)

#### REGION MIDDLE EAST, AFRICA

Drägerwerk AG & Co. KGaA  
Branch Office  
P.O. Box 505108  
Dubai, United Arab Emirates  
Tel +971 4 4294 600  
Fax +971 4 4294 699  
[contactuae@draeger.com](mailto:contactuae@draeger.com)

#### REGION ASIA / PACIFIC

Draeger Singapore Pte. Ltd.  
25 International Business Park  
#04-20/21 German Centre  
Singapore 609916  
Tel +65 6308 9400  
Fax +65 6308 9401  
[asia.pacific@draeger.com](mailto:asia.pacific@draeger.com)

#### REGION NORTH AMERICA

Draeger, Inc.  
3135 Quarry Road  
Telford, PA 18969-1042, USA  
Tel +1 800 4DRAGER  
(+1 800 437 2437)  
Fax +1 215 723 5935  
[info.usa@draeger.com](mailto:info.usa@draeger.com)

#### REGION CENTRAL AND SOUTH AMERICA

Dräger Panama Comercial  
S. de R.L.  
Complejo Business Park,  
V tower, 10th floor  
Panama City  
Tel +507 377 9100  
Fax +507 377 9130  
[contactcsa@draeger.com](mailto:contactcsa@draeger.com)