

OLYMPUS®

Your Vision, Our Future

TD-20

Manual Disinfector



The TD-20 Manual Disinfector, facilitates the cleaning and disinfection of Olympus flexible endoscopes.

Features:

- **Compact and mobile** – no requirement to be connected to water/ electricity supply.
- **Integral hand pump** – situated in the trolley leg for improved operation and aesthetics.
- **Three fluid containers with lids** – to allow scope processing from cleaning to disinfection and rinsing. The container lids minimise fluid evaporation, while the integral pouring lip and hand grip on the containers facilitate filling with little or no splashing.
- **Fluid irrigation connectors** – all connectors are integral to the top tray, minimising the risk of splashing.
- **Smoked polycarbonate lid** – helps prevent disinfectant evaporation while allowing visual monitoring of the scope at any stage of the cycle.
- **Integral timer with alarm** – allows the disinfectant contact time to be accurately monitored.
- **Instrument tray** – purpose designed to accept Olympus compatible endoscopes and their accessories in the minimum amount of disinfectant.


Construction:

The TD-20 Manual Disinfector has a powder-coated steel frame, vacuum-moulded ABS plastic instrument tray and shelf and smoked polycarbonate lid.

Three buckets with lids are supplied, together with channel connectors for compatible Olympus endoscopes and their accessories.

Occupational Health:

The Manual Disinfector requires the use of a cold fluid chemical disinfectant. It is important that whichever disinfectant is selected, the appropriate health and safety requirements for ventilation and personal protection are observed in accordance with the manufacturer's instructions.

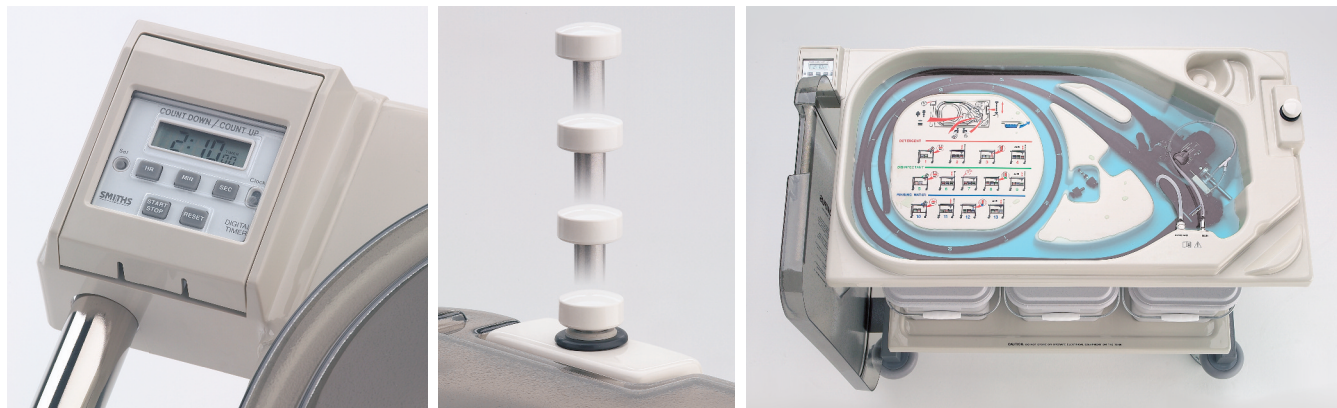

Specifications:

Dimensions: 960(W) x 930(H) x 510(D) mm

Weight: 31 kg

Instrument Tray Capacity: 6 litres

Olympus reserves the right to alter the specification without notice.

OLYMPUS[®]

OLYMPUS KEYMED GROUP
KeyMed House, Stock Road, Southend-on-Sea,
Essex SS2 5QH, United Kingdom
Telephone: +44 (0)1702 616333 Fax: +44 (0)1702 465677
e-mail: info@olympus.co.uk www.olympus.co.uk

OLYMPUS MEDICAL SYSTEMS CORP.
Shinjuku Monolith, 3-1, Nishi Shinjuku 2-chome, Shinjuku-ku, Tokyo 163-0914,

OLYMPUS EUROPA HOLDING GMBH
(Premises/Goods delivery) Wendenstraße 14-18, 20097 Hamburg, Germany
(Letters) Postfach 10 49 08, D-20034 Hamburg, Germany

OLYMPUS AMERICA INC.
3500 Corporate Parkway, PO Box 610, Center Valley, PA 18034-0610, USA

OLYMPUS SINGAPORE PTE LTD.
491B River Valley Road #12-01/04, Valley Point Office Tower, Singapore 248373, Singapore



ISO 9001 Registered
Certificate No. FM 20093
Quality Management System



ISO 14001 Registered
Certificate No. S465 00006
Environmental Management
System